

PROTEAN® Fragment Plating System

INSTRUCTIONS FOR USE (IFU) FOR THE UNITED STATES OF AMERICA

R: For use by physicians only. Federal Law restricts this device to sale by or on the order of a physician.

This package insert is designed to provide Instructions for Use of the PROTEAN® Fragment Plating System. **Failure to follow instructions may lead to patient injury.**

Description

The PROTEAN® Fragment Plating System is a set of titanium bone plates designed for stabilization and repair of small bone fragments. Included in the set are titanium and cobalt chrome bone screws, K-Wires, and specialized instrumentation. The system is provided non-sterile and is sterilized in the user facility.

The PROTEAN® Fragment Plating System is comprised of:

- Titanium alloy plates
- Titanium and cobalt chrome screws
- Stainless steel and cobalt chrome K-Wires (for provisional fixation not for implantation)

Indications

The PROTEAN® Fragment Plating System is intended for stabilization and fixation of small bone fragments in fresh fractures, revision procedures, joint fusion and reconstructions of small bones of the hand, foot, wrist, ankle, humerus, clavicle, scapula, finger, toe, and pelvis, particularly in osteopenic bone with the appropriate configuration of the device. The system is intended for use in adult and pediatric populations (children and adolescents).

Contraindications

Prior to using the PROTEAN® Fragment Plating System, ensure that none of the following patient conditions are present: active or latent infection, insufficient quantity or quality of bone and/or soft tissue, material sensitivity, or patients who are unwilling or incapable of following post operative care instructions. The use of the Protean Fragment Plating System is contraindicated in cases where the implant(s) would cross open epiphyseal plates in skeletally immature patients.

⚠ Warnings and Precautions

- If the screws are not attached and/or fully tightened, a non-union, delayed union or construct failure may occur.
- The information in this document should be shared with the patient.
- The patient should be informed about the importance of following the post operative rehabilitation prescribed to fully understand the limitations in activities of daily living. The patient must be warned that failure to follow postoperative care instructions may cause the implant or treatment to fail.
- Potential Fragment Plating System construct failures such as implant breakage, loosening, fixation failure, delayed union, or non-union may occur as a result of non-compliance to post operative rehabilitation, excessive wrist activities or construct overloading.
- DO NOT reuse any of the Fragment Plating System implantable components. Reuse may compromise the structural integrity of the construct.
DO NOT open the volar capsule as it may devascularize fracture fragments and destabilize the volar wrist ligaments.
- Protect the Fragment Plating System's implantable components against scratching or nicking. Such stress concentration can lead to implant failure.
- Before using the Fragment Plating System, inspect all implants and instruments for wear, disfiguration and physical damage. If evidence of wear, disfiguration or physical damage is found, DO NOT use and contact your local Skeletal Dynamics representative or the Skeletal Dynamics Customer Care Department.
- DO NOT permanently implant the Skeletal Dynamics K-Wires; they are intended to be used during provisional fixation of the Fragment Plating.
- DO NOT permanently implant the pre-loaded Drill Guides or AIMing Guides; they are intended to be removed prior to peg insertion.
- The Fragment Plating System has not been evaluated for safety and compatibility in the MR environment. It has not been tested for heating, migration, or image artifact in the MR environment. The safety of the Fragment

Plating System in the MR environment is unknown. Scanning a patient who has this device may result in patient injury.

- DO NOT use peg/screw lengths that will excessively protrude through the far cortex as it may result in soft tissue irritation.
- The maximum angulation of the PLS should not exceed 10° from the trajectory of the respective hole.
- The Non-locking Threaded Pegs are NOT intended to provide subchondral support. Their use should be limited to capture remote bone fragments where partially or fully threaded pegs cannot be used.
- The Fragment Plating System is to be used only with Skeletal Dynamics instruments, implants and accessories.
- Dispose of contaminated implants and instruments per established facility guidelines and protocols.
- Caution should be taken for interference to pacemakers during electrocautery or by uncertified drill power sources.
- Seek medical help immediately if implant malfunctions.
- To maintain traceability of the Fragment Plating System implantable components, you must record each of the respective components LOT numbers into the patient records.

Potential Adverse Events

The following are potential risks that have been associated with orthopedic surgery: infection, non-union, persistent pain, loosening or migration of the implants.

Directions for Use

The PROTEAN® Fragment Plating System should only be used by surgeons who have experience with this system. Each surgeon must evaluate the appropriateness for the use of the Fragment Plating System based on their clinical experiences. Please refer to the appropriate PROTEAN® Fragment Plating System's Surgical Technique Guide (MKT-00005-00, MKT-00082-00, MKT-00485-00 or MKT-00498-00) to review the surgical approach as described by Jorge L. Orbay, M.D. of the *Miami Hand and Upper Extremity Institute* located in Miami, Florida.

Cleaning, Sterilization, and Inspection

For instructions on cleaning, disinfection, sterilization, and inspection of the PROTEAN® Fragment Plating System products please refer to the cleaning and sterilization instructions for use (IFU-04056-00). Refer to the appropriate PROTEAN® Fragment Plating System Surgical Technique Guide (MKT-00005-00, MKT-00082-00, MKT-00485-00 or MKT-00498-00) for proper kitted tray arrangement.

Resources

The latest version of the Instructions for Use may be requested in a physical format by email (orders@skeletaldynamics.com) or by phone (+1-877-753-5396). The physical copy will be provided within 7 calendar days of receiving a request from the user or at the time of delivery of the device if so, requested at the time of order.

For the most current instructions for use visit www.skeletaldynamics.com/resources. Instructions for Use should always be reviewed before using or implanting a device.

Disclaimer of Warranty and Limited Remedies

Skeletal Dynamics, Inc makes no express or implied warranty, including any implied warranty of merchantability or fitness for a particular purpose, on the product(s) described in this publication. Skeletal Dynamics, Inc. shall not be liable under any circumstances for any direct, incidental or consequential damages other than as expressly provided by specific law. No person has authority to bind Skeletal Dynamics, Inc. to any representation or warranty except as specifically set forth in this publication. Descriptions or specifications provided by Skeletal Dynamics, Inc. in any publication is only included to generally describe the product when manufactured and do not constitute any express warranties.












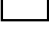























Skeletal Dynamics, Inc
7300 N. Kendall Dr. / Suite 800
Miami, FL 33156
United States of America
1-877-753-5396





















SYMBOLS GLOSSARY





Symbols that follow BS EN ISO 15223-1 / Medical Devices-Symbols to be used with medical device labels, labelling and information to be supplied.

Symbol	Symbol Reference Number and Title	Description
	5.2.7 Non-sterile	Indicates a medical device that has not been subjected to a sterilization process
	5.2.6 Do not re-sterilize	Indicates a medical device that is not to be re-sterilized
	5.1.6 Catalogue Number	Indicates the manufacturer's catalogue number so the medical device can be identified
	5.4.3 Consult Instructions for Use	Indicates the need for the user to consult the instructions for use
	5.1.5 Batch Code	Indicates the manufacturer's batch code or lot can be identified
	5.4.2 Do Not Re-use	Indicates a medical device that is intended for one single use only
	5.1.1 Manufacturer	Indicates the medical device manufacturer
	5.2.8 Do Not Use If Package is Damaged	Indicates a medical device that should not be used if the package has been damaged or opened
	5.1.4 Use-by-date	Indicates the date after which the medical device is not to be used
	5.7.10 Unique Device Identifier	Indicates a carrier that contains unique device identifier information
	5.7.7 Medical device	Indicated the item is a medical device
	5.1.3 Date of Manufacture	Indicates the date when the medical device was manufactured

PROTEAN Fragment Plate	
PROTEAN Fragment Plate, Double Hockey Stick PRT-FSP-LR (01)00841506102917  <small>(01) 00841506102917</small>	PROTEAN Fragment Plate, Y PRT-FSP-YS (01)00841506102931  <small>(01) 00841506102931</small>
PROTEAN Fragment Plate, Scapholunate, Short, Right PRT-FSP-SR1 (01)00841506103860  <small>(01) 00841506103860</small>	PROTEAN Fragment Plate, Scapholunate, Long, Right PRT-FSP-SR2 (01)00841506103884  <small>(01) 00841506103884</small>
PROTEAN Fragment Plate, Scapholunate, Short, Left PRT-FSP-SL1 (01)00841506103877  <small>(01) 00841506103877</small>	PROTEAN Fragment Plate, Scapholunate, Long, Left PRT-FSP-SL2 (01)00841506103891  <small>(01) 00841506103891</small>
PROTEAN Fragment Plate, Radial Column Plate, Right PRT-RCP-RT (01)00841506109923  <small>(01) 00841506109923</small>	PROTEAN Fragment Plate, Central Column Plate, Right PRT-CCP-RT (01)00841506109947  <small>(01) 00841506109947</small>
PROTEAN Fragment Plate, Radial Column Plate, Left PRT-RCP-LT (01)00841506109930  <small>(01) 00841506109930</small>	PROTEAN Fragment Plate, Central Column Plate, Left PRT-CCP-LT (01)00841506109954  <small>(01) 00841506109954</small>
PROTEAN Fragment Plate, Distal Ulna PRT-FSP-DU (01)00841506102900  <small>(01) 00841506102900</small>	
Threaded Peg, Locking Screws, 2.3mm	
Threaded Peg, Locking, 2.3mm x 10mm, Ti TPLS-23100-TS (01)00841506103358  <small>(01) 00841506103358</small>	Threaded Peg, Locking, 2.3mm x 20mm, Ti TPLS-23200-TS (01)00841506103426  <small>(01) 00841506103426</small>
Threaded Peg, Locking, 2.3mm x 12mm, Ti TPLS-23120-TS (01)00841506103365  <small>(01) 00841506103365</small>	Threaded Peg, Locking, 2.3mm x 21mm, Ti TPLS-23210-TS (01)00841506103433  <small>(01) 00841506103433</small>
Threaded Peg, Locking, 2.3mm x 14mm, Ti TPLS-23140-TS (01)00841506103372  <small>(01) 00841506103372</small>	Threaded Peg, Locking, 2.3mm x 22mm, Ti TPLS-23220-TS (01)00841506103440  <small>(01) 00841506103440</small>
Threaded Peg, Locking, 2.3mm x 16mm, Ti TPLS-23160-TS (01)00841506103389  <small>(01) 00841506103389</small>	Threaded Peg, Locking, 2.3mm x 23mm, Ti TPLS-23230-TS (01)00841506103457  <small>(01) 00841506103457</small>
Threaded Peg, Locking, 2.3mm x 17mm, Ti TPLS-23170-TS (01)00841506103396  <small>(01) 00841506103396</small>	Threaded Peg, Locking, 2.3mm x 24mm, Ti TPLS-23240-TS (01)00841506103464  <small>(01) 00841506103464</small>





















Threaded Peg, Locking, 2.3mm x 18mm, Ti TPLS-23180-TS (01)00841506103402  <small>(01) 00841506103402</small>	Threaded Peg, Locking, 2.3mm x 26mm, Ti TPLS-23260-TS (01)00841506103471  <small>(01) 00841506103471</small>
Threaded Peg, Locking, 2.3mm x 19mm, Ti TPLS-23190-TS (01)00841506103419  <small>(01) 00841506103419</small>	Threaded Peg, Locking, 2.3mm x 28mm, Ti TPLS-23280-TS (01)00841506103488  <small>(01) 00841506103488</small>
High Compression Locking Screws, 2.7mm	
High Compression Locking Peg, 2.7mm x 10mm, Ti HCLP-27100-TS (01)00841506101682  <small>(01) 00841506101682</small>	High Compression Locking Peg, 2.7mm x 21mm, Ti HCLP-27210-TS (01)00841506101750  <small>(01) 00841506101750</small>
High Compression Locking Peg, 2.7mm x 12mm, Ti HCLP-27120-TS (01)00841506101699  <small>(01) 00841506101699</small>	High Compression Locking Peg, 2.7mm x 22mm, Ti HCLP-27220-TS (01)00841506101767  <small>(01) 00841506101767</small>
High Compression Locking Peg, 2.7mm x 14mm, Ti HCLP-27140-TS (01)00841506101705  <small>(01) 00841506101705</small>	High Compression Locking Peg, 2.7mm x 23mm, Ti HCLP-27230-TS (01)00841506101774  <small>(01) 00841506101774</small>
High Compression Locking Peg, 2.7mm x 16mm, Ti HCLP-27160-TS (01)00841506101712  <small>(01) 00841506101712</small>	High Compression Locking Peg, 2.7mm x 24mm, Ti HCLP-27240-TS (01)00841506101781  <small>(01) 00841506101781</small>
High Compression Locking Peg, 2.7mm x 18mm, Ti HCLP-27180-TS (01)00841506101729  <small>(01) 00841506101729</small>	High Compression Locking Peg, 2.7mm x 26mm, Ti HCLP-27260-TS (01)00841506101798  <small>(01) 00841506101798</small>
High Compression Locking Peg, 2.7mm x 19mm, Ti HCLP-27190-TS (01)00841506101736  <small>(01) 00841506101736</small>	High Compression Locking Peg, 2.7mm x 28mm, Ti HCLP-27280-TS (01)00841506101804  <small>(01) 00841506101804</small>
High Compression Locking Peg, 2.7mm x 20mm, Ti HCLP-27200-TS (01)00841506101743  <small>(01) 00841506101743</small>	High Compression Locking Peg, 2.7mm x 30mm, Ti HCLP-27300-TS (01)00841506101811  <small>(01) 00841506101811</small>
Threaded Pegs, Non-Locking, 2.7mm	
Threaded Peg, Non-Locking, 2.7mm x 10mm, Ti TPNL-27100-TS (01)00841506103518  <small>(01) 00841506103518</small>	Threaded Peg, Non-Locking, 2.7mm x 22mm, Ti TPNL-27220-TS (01)00841506103570  <small>(01) 00841506103570</small>
Threaded Peg, Non-Locking, 2.7mm x 12mm, Ti TPNL-27120-TS (01)00841506103525  <small>(01) 00841506103525</small>	Threaded Peg, Non-Locking, 2.7mm x 24mm, Ti TPNL-27240-TS (01)00841506103587  <small>(01) 00841506103587</small>






















Threaded Peg, Non-Locking, 2.7mm x 14mm, Ti TPNL-27140-TS (01)00841506103532  <small>(01) 00841506103532</small>	Threaded Peg, Non-Locking, 2.7mm x 26mm, Ti TPNL-27260-TS (01)00841506103594  <small>(01) 00841506103594</small>
Threaded Peg, Non-Locking, 2.7mm x 16mm, Ti TPNL-27160-TS (01)00841506103549  <small>(01) 00841506103549</small>	Threaded Peg, Non-Locking, 2.7mm x 28mm, Ti TPNL-27280-TS (01)00841506103600  <small>(01) 00841506103600</small>
Threaded Peg, Non-Locking, 2.7mm x 18mm, Ti TPNL-27180-TS (01)00841506103556  <small>(01) 00841506103556</small>	Threaded Peg, Non-Locking, 2.7mm x 30mm, Ti TPNL-27300-TS (01)00841506103617  <small>(01) 00841506103617</small>
Threaded Peg, Non-Locking, 2.7mm x 20mm, Ti TPNL-27200-TS (01)00841506103563  <small>(01) 00841506103563</small>	
Polyaxial Locking Screws, 2.5mm	
Screw, Polyaxial Locking, 2.5mm x 10mm, CoCr PALS-25100-CC (01)00841506102665  <small>(01) 00841506102665</small>	Screw, Polyaxial Locking, 2.5mm x 20mm, CoCr PALS-25200-CC (01)00841506102719  <small>(01) 00841506102719</small>
Screw, Polyaxial Locking, 2.5mm x 12mm, CoCr PALS-25120-CC (01)00841506102672  <small>(01) 00841506102672</small>	Screw, Polyaxial Locking, 2.5mm x 22mm, CoCr PALS-25220-CC (01)00841506102726  <small>(01) 00841506102726</small>
Screw, Polyaxial Locking, 2.5mm x 14mm, CoCr PALS-25140-CC (01)00841506102689  <small>(01) 00841506102689</small>	Screw, Polyaxial Locking, 2.5mm x 24mm, CoCr PALS-25240-CC (01)00841506102733  <small>(01) 00841506102733</small>
Screw, Polyaxial Locking, 2.5mm x 16mm, CoCr PALS-25160-CC (01)00841506102696  <small>(01) 00841506102696</small>	Screw, Polyaxial Locking, 2.5mm x 26mm, CoCr PALS-25260-CC (01)00841506102740  <small>(01) 00841506102740</small>
Screw, Polyaxial Locking, 2.5mm x 18mm, CoCr PALS-25180-CC (01)00841506102702  <small>(01) 00841506102702</small>	Screw, Polyaxial Locking, 2.5mm x 28mm, CoCr PALS-25280-CC (01)00841506102757  <small>(01) 00841506102757</small>
Washer Button	
Washer, Button WBTN-2750-T (01)00841506103730  <small>(01) 00841506103730</small>	

Guides			
AIMing Guides, 1.5mm PDG-AIM-015 (01)00841506102870	 <small>(01) 00841506102870</small>	PLS AIMing Guides, .9mm x 10° PLS-AIM-0910 (01)00841506102887	 <small>(01) 00841506102887</small>
K-Wires			
K-Wire, Standard Tip, .9mm x 152mm KWIR-STD-09152 (01)00841506102498	 <small>(01) 00841506102498</small>	K-Wire, 1.6mm x 127mm KWIR-STD-15127 (01)00841506102504	 <small>(01) 00841506102504</small>

K-Wires	
K-Wire, Single Diamond, 0.9mm x 89mm, CoCr KWCC-SD-09089 (01)00841506137544  <small>(01) 00841506137544</small>	K-Wire, Single Diamond, 1.2mm x 89mm CoCr KWCC-SD-12089 (01)00841506138107  <small>(01) 00841506138107</small>
K-Wire, Double Diamond, 0.9mm x 89mm, CoCr KWCC-DT-09089 (01)00841506137551  <small>(01) 00841506137551</small>	K-Wire, Double Double, 1.2mm x 89mm CoCr KWCC-DD-12089 (01)00841506138114  <small>(01) 00841506138114</small>
K-Wire, 1.6mm x 127mm KWIR-STD-15127 (01)00841506102504  <small>(01) 00841506102504</small>	K-Wire, Threaded Tip with Olive Stop, 1.2mm x 64mm KWIR-OLV-12064 (01)00841506146690  <small>(01) 00841506146690</small>
AIMing Guides	
PROTEAN AIMing Guide, 2.0mm PRT-AIM-20 (01)00841506138091  <small>(01) 00841506138091</small>	PROTEAN AIMing Guide, 2.5mm PRT-AIM-25 (01)00841506138381  <small>(01) 00841506138381</small>
Non-Locking Screws, 1.2mm	
Screw, Non-Locking, 1.2mm x 6mm, CoCr VAN-12006-CC (01)00841506137216  <small>(01) 00841506137216</small>	Screw, Non-Locking, 1.2mm x 11mm, CoCr VAN-12011-CC (01)00841506137261  <small>(01) 00841506137261</small>
Screw, Non-Locking, 1.2mm x 7mm, CoCr VAN-12007-CC (01)00841506137223  <small>(01) 00841506137223</small>	Screw, Non-Locking, 1.2mm x 12mm, CoCr VAN-12012-CC (01)00841506137278  <small>(01) 00841506137278</small>
Screw, Non-Locking, 1.2mm x 8mm, CoCr VAN-12008-CC (01)00841506137230  <small>(01) 00841506137230</small>	Screw, Non-Locking, 1.2mm x 13mm, CoCr VAN-12013-CC (01)00841506137285  <small>(01) 00841506137285</small>
Screw, Non-Locking, 1.2mm x 9mm, CoCr VAN-12009-CC (01)00841506137247  <small>(01) 00841506137247</small>	Screw, Non-Locking, 1.2mm x 14mm, CoCr VAN-12014-CC (01)00841506137292  <small>(01) 00841506137292</small>
Screw, Non-Locking, 1.2mm x 10mm, CoCr VAN-12010-CC (01)00841506137254  <small>(01) 00841506137254</small>	
PROTEAN Hand Plates, 1.5mm	
PROTEAN Fragment Plate, Ladder, 1.5mm, Ti PRT-LDR-015 (01)00841506137193  <small>(01) 00841506137193</small>	PROTEAN Fragment Plate, Straight/Rotation, 1.5mm, Ti PRT-STR-015 (01)00841506137186  <small>(01) 00841506137186</small>

PROTEAN Fragment Plate, Phalanx, 1.5, Ti PRT-PHX-015 (01)00841506137179  <small>(01) 00841506137179</small>	PROTEAN Fragment Plate, Web, 1.5, Ti PRT-WEB-015 (01)00841506137209  <small>(01) 00841506137209</small>
PROTEAN Fragment Plate, TY, 1.5, Ti PRT-TY-015 (01)00841506137162  <small>(01) 00841506137162</small>	PROTEAN Fragment Plate, Metacarpal, 1.5, Ti PRT-MC-015 (01)00841506145228  <small>(01) 00841506145228</small>
Variable Angle Locking Screws, 1.5mm	
Screw, Variable Angle Locking, 1.5mm x 5mm, CoCr VAL-15050-CC (01)00841506137315  <small>(01) 00841506137315</small>	Screw, Variable Angle Locking, 1.5mm x 15mm, CoCr VAL-15150-CC (01)00841506136479  <small>(01) 00841506136479</small>
Screw, Variable Angle Locking, 1.5mm x 6mm, CoCr VAL-15060-CC (01)00841506136387  <small>(01) 00841506136387</small>	Screw, Variable Angle Locking, 1.5mm x 16mm, CoCr VAL-15160-CC (01)00841506136486  <small>(01) 00841506136486</small>
Screw, Variable Angle Locking, 1.5mm x 7mm, CoCr VAL-15070-CC (01)00841506136394  <small>(01) 00841506136394</small>	Screw, Variable Angle Locking, 1.5mm x 18mm, CoCr VAL-15180-CC (01)00841506136493  <small>(01) 00841506136493</small>
Screw, Variable Angle Locking, 1.5mm x 8mm, CoCr VAL-15080-CC (01)00841506136400  <small>(01) 00841506136400</small>	Screw, Variable Angle Locking, 1.5mm x 20mm, CoCr VAL-15200-CC (01)00841506136509  <small>(01) 00841506136509</small>
Screw, Variable Angle Locking, 1.5mm x 9mm, CoCr VAL-15090-CC (01)00841506136417  <small>(01) 00841506136417</small>	Screw, Variable Angle Locking, 1.5mm x 22mm, CoCr VAL-15220-CC (01)00841506136516  <small>(01) 00841506136516</small>
Screw, Variable Angle Locking, 1.5mm x 10mm, CoCr VAL-15100-CC (01)00841506136424  <small>(01) 00841506136424</small>	Screw, Variable Angle Locking, 1.5mm x 24mm, CoCr VAL-15240-CC (01)00841506136523  <small>(01) 00841506136523</small>
Screw, Variable Angle Locking, 1.5mm x 11mm, CoCr VAL-15110-CC (01)00841506136431  <small>(01) 00841506136431</small>	Screw, Variable Angle Locking, 1.5mm x 26mm, CoCr VAL-15260-CC (01)00841506136530  <small>(01) 00841506136530</small>
Screw, Variable Angle Locking, 1.5mm x 12mm, CoCr VAL-15120-CC (01)00841506136448  <small>(01) 00841506136448</small>	Screw, Variable Angle Locking, 1.5mm x 28mm, CoCr VAL-15280-CC (01)00841506136547  <small>(01) 00841506136547</small>
Screw, Variable Angle Locking, 1.5mm x 13mm, CoCr VAL-15130-CC (01)00841506136455  <small>(01) 00841506136455</small>	Screw, Variable Angle Locking, 1.5mm x 30mm, CoCr VAL-15300-CC (01)00841506136554  <small>(01) 00841506136554</small>

Screw, Variable Angle Locking, 1.5mm x 14mm, CoCr VAL-15140-CC (01)00841506136462  <small>(01) 00841506136462</small>	
Locking Screws, 1.5mm	
Screw, Locking, 1.5mm x 5mm, Ti FAL-15050-TS (01)00841506146904  <small>(01) 00841506146904</small>	Screw, Locking, 1.5mm x 15mm, Ti FAL-15150-TS (01)00841506147000  <small>(01) 00841506147000</small>
Screw, Locking, 1.5mm x 6mm, Ti FAL-15060-TS (01)00841506146911  <small>(01) 00841506146911</small>	Screw, Locking, 1.5mm x 16mm, Ti FAL-15160-TS (01)00841506147017  <small>(01) 00841506147017</small>
Screw, Locking, 1.5mm x 7mm, Ti FAL-15070-TS (01)00841506146928  <small>(01) 00841506146928</small>	Screw, Locking, 1.5mm x 18mm, Ti FAL-15180-TS (01)00841506147024  <small>(01) 00841506147024</small>
Screw, Locking, 1.5mm x 8mm, Ti FAL-15080-TS (01)00841506146935  <small>(01) 00841506146935</small>	Screw, Locking, 1.5mm x 20mm, Ti FAL-15200-TS (01)00841506147031  <small>(01) 00841506147031</small>
Screw, Locking, 1.5mm x 9mm, Ti FAL-15090-TS (01)00841506146942  <small>(01) 00841506146942</small>	Screw, Locking, 1.5mm x 22mm, Ti FAL-15220-TS (01)00841506147048  <small>(01) 00841506147048</small>
Screw, Locking, 1.5mm x 10mm, Ti FAL-15100-TS (01)00841506146959  <small>(01) 00841506146959</small>	Screw, Locking, 1.5mm x 24mm, Ti FAL-15240-TS (01)00841506147055  <small>(01) 00841506147055</small>
Screw, Locking, 1.5mm x 11mm, Ti FAL-15110-TS (01)00841506146966  <small>(01) 00841506146966</small>	Screw, Locking, 1.5mm x 26mm, Ti FAL-15260-TS (01)00841506147062  <small>(01) 00841506147062</small>
Screw, Locking, 1.5mm x 12mm, Ti FAL-15120-TS (01)00841506146973  <small>(01) 00841506146973</small>	Screw, Locking, 1.5mm x 28mm, Ti FAL-15280-TS (01)00841506147079  <small>(01) 00841506147079</small>
Screw, Locking, 1.5mm x 13mm, Ti FAL-15130-TS (01)00841506146980  <small>(01) 00841506146980</small>	Screw, Locking, 1.5mm x 30mm, Ti FAL-15300-TS (01)00841506147086  <small>(01) 00841506147086</small>
Screw, Locking, 1.5mm x 14mm, Ti FAL-15140-TS (01)00841506146997  <small>(01) 00841506146997</small>	

Non-Locking Screws, 1.5mm	
Screw, Non-Locking, 1.5mm x 5mm, Ti VAN-15050-TS (01)00841506137339  <small>(01) 00841506137339</small>	Screw, Non-Locking, 1.5mm x 15mm, Ti VAN-15150-TS (01)00841506137438  <small>(01) 00841506137438</small>
Screw, Non-Locking, 1.5mm x 6mm, Ti VAN-15060-TS (01)00841506137346  <small>(01) 00841506137346</small>	Screw, Non-Locking, 1.5mm x 16mm, Ti VAN-15160-TS (01)00841506137445  <small>(01) 00841506137445</small>
Screw, Non-Locking, 1.5mm x 7mm, Ti VAN-15070-TS (01)00841506137353  <small>(01) 00841506137353</small>	Screw, Non-Locking, 1.5mm x 18mm, Ti VAN-15180-TS (01)00841506137452  <small>(01) 00841506137452</small>
Screw, Non-Locking, 1.5mm x 8mm, Ti VAN-15080-TS (01)00841506137360  <small>(01) 00841506137360</small>	Screw, Non-Locking, 1.5mm x 20mm, Ti VAN-15200-TS (01)00841506137469  <small>(01) 00841506137469</small>
Screw, Non-Locking, 1.5mm x 9mm, Ti VAN-15090-TS (01)00841506137377  <small>(01) 00841506137377</small>	Screw, Non-Locking, 1.5mm x 22mm, Ti VAN-15220-TS (01)00841506137476  <small>(01) 00841506137476</small>
Screw, Non-Locking, 1.5mm x 10mm, Ti VAN-15100-TS (01)00841506137384  <small>(01) 00841506137384</small>	Screw, Non-Locking, 1.5mm x 24mm, Ti VAN-15240-TS (01)00841506137483  <small>(01) 00841506137483</small>
Screw, Non-Locking, 1.5mm x 11mm, Ti VAN-15110-TS (01)00841506137391  <small>(01) 00841506137391</small>	Screw, Non-Locking, 1.5mm x 26mm, Ti VAN-15260-TS (01)00841506137490  <small>(01) 00841506137490</small>
Screw, Non-Locking, 1.5mm x 12mm, Ti VAN-15120-TS (01)00841506137407  <small>(01) 00841506137407</small>	Screw, Non-Locking, 1.5mm x 28mm, Ti VAN-15280-TS (01)00841506137506  <small>(01) 00841506137506</small>
Screw, Non-Locking, 1.5mm x 13mm, Ti VAN-15130-TS (01)00841506137414  <small>(01) 00841506137414</small>	Screw, Non-Locking, 1.5mm x 30mm, Ti VAN-15300-TS (01)00841506137513  <small>(01) 00841506137513</small>
Screw, Non-Locking, 1.5mm x 14mm, Ti VAN-15140-TS (01)00841506137421  <small>(01) 00841506137421</small>	
PROTEAN Hand Plates, 2.0mm	
PROTEAN Fragment Plate, TY, 2.0mm, Ti PRT-TY-020 (01)00841506137605  <small>(01) 00841506137605</small>	PROTEAN Fragment Plate, Straight/Rotation, 2.0mm, Ti PRT-STR-020 (01)00841506137636  <small>(01) 00841506137636</small>

PROTEAN Fragment Plate, Metacarpal, 2.0mm, Ti PRT-MC-020 (01)00841506137612  <small>(01) 00841506137612</small>	PROTEAN Fragment Plate, Ladder, 2.0mm, Ti PRT-LDR-020 (01)00841506137643  <small>(01) 00841506137643</small>
PROTEAN Fragment Plate, Phalanx, 2.0mm, Ti PRT-PHX-020 (01)00841506137629  <small>(01) 00841506137629</small>	PROTEAN Fragment Plate, Web, 2.0mm, Ti PRT-WEB-020 (01)00841506137650  <small>(01) 00841506137650</small>
PROTEAN Fragment Plate, PS, 2.0, Ti PRT-PS-020 (01)00841506140346  <small>(01) 00841506140346</small>	PROTEAN Fragment Plate, PS, Large, 2.0, Ti PRT-PSO-020 (01)00841506140353  <small>(01) 00841506140353</small>
Variable Angle Locking Screws, 2.0mm	
Screw, Variable Angle Locking, 2.0mm x 6mm, CoCr VAL-20060-CC (01)00841506136202  <small>(01) 00841506136202</small>	Screw, Variable Angle Locking, 2.0mm x 15mm, CoCr VAL-20150-CC (01)00841506136295  <small>(01) 00841506136295</small>
Screw, Variable Angle Locking, 2.0mm x 7mm, CoCr VAL-20070-CC (01)00841506136219  <small>(01) 00841506136219</small>	Screw, Variable Angle Locking, 2.0mm x 16mm, CoCr VAL-20160-CC (01)00841506136301  <small>(01) 00841506136301</small>
Screw, Variable Angle Locking, 2.0mm x 8mm, CoCr VAL-20080-CC (01)00841506136226  <small>(01) 00841506136226</small>	Screw, Variable Angle Locking, 2.0mm x 18mm, CoCr VAL-20180-CC (01)00841506136318  <small>(01) 00841506136318</small>
Screw, Variable Angle Locking, 2.0mm x 9mm, CoCr VAL-20090-CC (01)00841506136233  <small>(01) 00841506136233</small>	Screw, Variable Angle Locking, 2.0mm x 20mm, CoCr VAL-20200-CC (01)00841506136325  <small>(01) 00841506136325</small>
Screw, Variable Angle Locking, 2.0mm x 10mm, CoCr VAL-20100-CC (01)00841506136240  <small>(01) 00841506136240</small>	Screw, Variable Angle Locking, 2.0mm x 22mm, CoCr VAL-20220-CC (01)00841506136332  <small>(01) 00841506136332</small>
Screw, Variable Angle Locking, 2.0mm x 11mm, CoCr VAL-20110-CC (01)00841506136257  <small>(01) 00841506136257</small>	Screw, Variable Angle Locking, 2.0mm x 24mm, CoCr VAL-20240-CC (01)00841506136349  <small>(01) 00841506136349</small>
Screw, Variable Angle Locking, 2.0mm x 12mm, CoCr VAL-20120-CC (01)00841506136264  <small>(01) 00841506136264</small>	Screw, Variable Angle Locking, 2.0mm x 26mm, CoCr VAL-20260-CC (01)00841506136356  <small>(01) 00841506136356</small>
Screw, Variable Angle Locking, 2.0mm x 13mm, CoCr VAL-20130-CC (01)00841506136271  <small>(01) 00841506136271</small>	Screw, Variable Angle Locking, 2.0mm x 28mm, CoCr VAL-20280-CC (01)00841506136363  <small>(01) 00841506136363</small>

Screw, Variable Angle Locking, 2.0mm x 14mm, CoCr VAL-20140-CC (01)00841506136288  <small>(01) 00841506136288</small>	Screw, Variable Angle Locking, 2.0mm x 30mm, CoCr VAL-20300-CC (01)00841506136370  <small>(01) 00841506136370</small>
Locking Screws, 2.0mm	
Screw, Locking, 2.0mm x 6mm, Ti FAL-20060-TS (01)00841506147093  <small>(01) 00841506147093</small>	Screw, Locking, 2.0mm x 15mm, Ti FAL-20150-TS (01)00841506147215  <small>(01) 00841506147215</small>
Screw, Locking, 2.0mm x 7mm, Ti FAL-20070-TS (01)00841506147109  <small>(01) 00841506147109</small>	Screw, Locking, 2.0mm x 16mm, Ti FAL-20160-TS (01)00841506147192  <small>(01) 00841506147192</small>
Screw, Locking, 2.0mm x 8mm, Ti FAL-20080-TS (01)00841506147116  <small>(01) 00841506147116</small>	Screw, Locking, 2.0mm x 18mm, Ti FAL-20180-TS (01)00841506147208  <small>(01) 00841506147208</small>
Screw, Locking, 2.0mm x 9mm, Ti FAL-20090-TS (01)00841506147123  <small>(01) 00841506147123</small>	Screw, Locking, 2.0mm x 20mm, Ti FAL-20200-TS (01)00841506147215  <small>(01) 00841506147215</small>
Screw, Locking, 2.0mm x 10mm, Ti FAL-20100-TS (01)00841506147130  <small>(01) 00841506147130</small>	Screw, Locking, 2.0mm x 22mm, Ti FAL-20220-TS (01)00841506147222  <small>(01) 00841506147222</small>
Screw, Locking, 2.0mm x 11mm, Ti FAL-20110-TS (01)00841506147147  <small>(01) 00841506147147</small>	Screw, Locking, 2.0mm x 24mm, Ti FAL-20240-TS (01)00841506147239  <small>(01) 00841506147239</small>
Screw, Locking, 2.0mm x 12mm, Ti FAL-20120-TS (01)00841506147154  <small>(01) 00841506147154</small>	Screw, Locking, 2.0mm x 26mm, Ti FAL-20260-TS (01)00841506147246  <small>(01) 00841506147246</small>
Screw, Locking, 2.0mm x 13mm, Ti FAL-20130-TS (01)00841506147161  <small>(01) 00841506147161</small>	Screw, Locking, 2.0mm x 28mm, Ti FAL-20280-TS (01)00841506147253  <small>(01) 00841506147253</small>
Screw, Locking, 2.0mm x 14mm, Ti FAL-20140-TS (01)00841506147178  <small>(01) 00841506147178</small>	Screw, Locking, 2.0mm x 30mm, Ti FAL-20300-TS (01)00841506147260  <small>(01) 00841506147260</small>
Non-Locking Screws, 2.0mm	
Screw, Non-Locking, 2.0mm x 6mm, Ti VAN-20060-TS (01)00841506137667  <small>(01) 00841506137667</small>	Screw, Non-Locking, 2.0mm x 15mm, Ti VAN-20150-TS (01)00841506137759  <small>(01) 00841506137759</small>

Screw, Non-Locking, 2.0mm x 7mm, Ti VAN-20070-TS (01)00841506137674  <small>(01) 00841506137674</small>	Screw, Non-Locking, 2.0mm x 16mm, Ti VAN-20160-TS (01)00841506137766  <small>(01) 00841506137766</small>
Screw, Non-Locking, 2.0mm x 8mm, Ti VAN-20080-TS (01)00841506137681  <small>(01) 00841506137681</small>	Screw, Non-Locking, 2.0mm x 18mm, Ti VAN-20180-TS (01)00841506137773  <small>(01) 00841506137773</small>
Screw, Non-Locking, 2.0mm x 9mm, Ti VAN-20090-TS (01)00841506137698  <small>(01) 00841506137698</small>	Screw, Non-Locking, 2.0mm x 20mm, Ti VAN-20200-TS (01)00841506137780  <small>(01) 00841506137780</small>
Screw, Non-Locking, 2.0mm x 10mm, Ti VAN-20100-TS (01)00841506137704  <small>(01) 00841506137704</small>	Screw, Non-Locking, 2.0mm x 22mm, Ti VAN-20220-TS (01)00841506137797  <small>(01) 00841506137797</small>
Screw, Non-Locking, 2.0mm x 11mm, Ti VAN-20110-TS (01)00841506137711  <small>(01) 00841506137711</small>	Screw, Non-Locking, 2.0mm x 24mm, Ti VAN-20240-TS (01)00841506137803  <small>(01) 00841506137803</small>
Screw, Non-Locking, 2.0mm x 12mm, Ti VAN-20120-TS (01)00841506137728  <small>(01) 00841506137728</small>	Screw, Non-Locking, 2.0mm x 26mm, Ti VAN-20260-TS (01)00841506137810  <small>(01) 00841506137810</small>
Screw, Non-Locking, 2.0mm x 13mm, Ti VAN-20130-TS (01)00841506137735  <small>(01) 00841506137735</small>	Screw, Non-Locking, 2.0mm x 28mm, Ti VAN-20280-TS (01)00841506137827  <small>(01) 00841506137827</small>
Screw, Non-Locking, 2.0mm x 14mm, Ti VAN-20140-TS (01)00841506137742  <small>(01) 00841506137742</small>	Screw, Non-Locking, 2.0mm x 30mm, Ti VAN-20300-TS (01)00841506137834  <small>(01) 00841506137834</small>
Protean Hand Plates, 2.5mm	
PROTEAN Fragment Plate, TY, 2.5mm, Ti PRT-TY-025 (01)00841506137865  <small>(01) 00841506137865</small>	PROTEAN Fragment Plate, Straight/Rotation, 2.5mm, Ti PRT-STR-025 (01)00841506137889  <small>(01) 00841506137889</small>
PROTEAN Fragment Plate, Metacarpal, 2.5mm, Ti PRT-MC-025 (01)00841506137872  <small>(01) 00841506137872</small>	PROTEAN Fragment Plate, CMC PRT-MFP-CMC (01)00841506113050  <small>(01) 00841506113050</small>
Variable Angle Locking Screws, 2.5mm	
Screw, Variable Angle Locking, 2.5mm x 6mm, CoCr VAL-25060-CC (01)00841506136028  <small>(01) 00841506136028</small>	Screw, Variable Angle Locking, 2.5mm x 15mm, CoCr VAL-25150-CC (01)00841506136110  <small>(01) 00841506136110</small>

Screw, Variable Angle Locking, 2.5mm x 7mm, CoCr VAL-25070-CC (01)00841506136035  <small>(01) 00841506136035</small>	Screw, Variable Angle Locking, 2.5mm x 16mm, CoCr VAL-25160-CC (01)00841506136127  <small>(01) 00841506136127</small>
Screw, Variable Angle Locking, 2.5mm x 8mm, CoCr VAL-25080-CC (01)00841506136042  <small>(01) 00841506136042</small>	Screw, Variable Angle Locking, 2.5mm x 18mm, CoCr VAL-25180-CC (01)00841506136134  <small>(01) 00841506136134</small>
Screw, Variable Angle Locking, 2.5mm x 9mm, CoCr VAL-25090-CC (01)00841506136059  <small>(01) 00841506136059</small>	Screw, Variable Angle Locking, 2.5mm x 20mm, CoCr VAL-25200-CC (01)00841506136141  <small>(01) 00841506136141</small>
Screw, Variable Angle Locking, 2.5mm x 10mm, CoCr VAL-25100-CC (01)00841506136066  <small>(01) 00841506136066</small>	Screw, Variable Angle Locking, 2.5mm x 22mm, CoCr VAL-25220-CC (01)00841506136158  <small>(01) 00841506136158</small>
Screw, Variable Angle Locking, 2.5mm x 11mm, CoCr VAL-25110-CC (01)00841506136073  <small>(01) 00841506136073</small>	Screw, Variable Angle Locking, 2.5mm x 24mm, CoCr VAL-25240-CC (01)00841506136165  <small>(01) 00841506136165</small>
Screw, Variable Angle Locking, 2.5mm x 12mm, CoCr VAL-25120-CC (01)00841506136080  <small>(01) 00841506136080</small>	Screw, Variable Angle Locking, 2.5mm x 26mm, CoCr VAL-25260-CC (01)00841506136172  <small>(01) 00841506136172</small>
Screw, Variable Angle Locking, 2.5mm x 13mm, CoCr VAL-25130-CC (01)00841506136097  <small>(01) 00841506136097</small>	Screw, Variable Angle Locking, 2.5mm x 28mm, CoCr VAL-25280-CC (01)00841506136189  <small>(01) 00841506136189</small>
Screw, Variable Angle Locking, 2.5mm x 14mm, CoCr VAL-25140-CC (01)00841506136103  <small>(01) 00841506136103</small>	Screw, Variable Angle Locking, 2.5mm x 30mm, CoCr VAL-25300-CC (01)00841506136196  <small>(01) 00841506136196</small>
Locking Screw, 2.5mm	
Screw, Locking, 2.5mm x 6mm, Ti FAL-25060-TS (01)00841506147277  <small>(01) 00841506147277</small>	Screw, Locking, 2.5mm x 15mm, Ti FAL-25150-TS (01)00841506147369  <small>(01) 00841506147369</small>
Screw, Locking, 2.5mm x 7mm, Ti FAL-25070-TS (01)00841506147284  <small>(01) 00841506147284</small>	Screw, Locking, 2.5mm x 16mm, Ti FAL-25160-TS (01)00841506147376  <small>(01) 00841506147376</small>
Screw, Locking, 2.5mm x 8mm, Ti FAL-25080-TS (01)00841506147291  <small>(01) 00841506147291</small>	Screw, Locking, 2.5mm x 18mm, Ti FAL-25180-TS (01)00841506147383  <small>(01) 00841506147383</small>
Screw, Locking, 2.5mm x 9mm, Ti FAL-25090-TS (01)00841506147307  <small>(01) 00841506147307</small>	IScrew, Locking, 2.5mm x 20mm, Ti FAL-25200-TS (01)00841506147390  <small>(01) 00841506147390</small>

Screw, Locking, 2.5mm x 10mm, Ti FAL-25100-TS (01)00841506147314	 (01) 00841506147314	Screw, Locking, 2.5mm x 22mm, Ti FAL-25220-TS (01)00841506147406	 (01) 00841506147406
Screw, Locking, 2.5mm x 11mm, Ti FAL-25110-TS (01)00841506147321	 (01) 00841506147321	Screw, Locking, 2.5mm x 24mm, Ti FAL-25240-TS (01)00841506147413	 (01) 00841506147413
Screw, Locking, 2.5mm x 12mm, Ti FAL-25120-TS (01)00841506147338	 (01) 00841506147338	Screw, Locking, 2.5mm x 26mm, Ti FAL-25260-TS (01)00841506147420	 (01) 00841506147420
Screw, Locking, 2.5mm x 13mm, Ti FAL-25130-TS (01)00841506147345	 (01) 00841506147345	Screw, Locking, 2.5mm x 28mm, Ti FAL-25280-TS (01)00841506147437	 (01) 00841506147437
Screw, Locking, 2.5mm x 14mm, Ti FAL-25140-TS (01)00841506147352	 (01) 00841506147352	Screw, Locking, 2.5mm x 30mm, Ti FAL-25300-TS (01)00841506147444	 (01) 00841506147444
Non-Locking Screw, 2.5mm			
Screw, Non-Locking, 2.5mm x 6mm, Ti VAN-25060-TS (01)00841506137896	 (01) 00841506137896	Screw, Non-Locking, 2.5mm x 15mm, Ti VAN-25150-TS (01)00841506137988	 (01) 00841506137988
Screw, Non-Locking, 2.5mm x 7mm, Ti VAN-25070-TS (01)00841506137902	 (01) 00841506137902	Screw, Non-Locking, 2.5mm x 16mm, Ti VAN-25160-TS (01)00841506137995	 (01) 00841506137995
Screw, Non-Locking, 2.5mm x 8mm, Ti VAN-25080-TS (01)00841506137919	 (01) 00841506137919	Screw, Non-Locking, 2.5mm x 18mm, Ti VAN-25180-TS (01)00841506138008	 (01) 00841506138008
Screw, Non-Locking, 2.5mm x 9mm, Ti VAN-25090-TS (01)00841506137926	 (01) 00841506137926	Screw, Non-Locking, 2.5mm x 20mm, Ti VAN-25200-TS (01)00841506138015	 (01) 00841506138015
Screw, Non-Locking, 2.5mm x 10mm, Ti VAN-25100-TS (01)00841506137933	 (01) 00841506137933	Screw, Non-Locking, 2.5mm x 22mm, Ti VAN-25220-TS (01)00841506138022	 (01) 00841506138022
Screw, Non-Locking, 2.5mm x 11mm, Ti VAN-25110-TS (01)00841506137940	 (01) 00841506137940	Screw, Non-Locking, 2.5mm x 24mm, Ti VAN-25240-TS (01)00841506138039	 (01) 00841506138039
Screw, Non-Locking, 2.5mm x 12mm, Ti VAN-25120-TS (01)00841506137957	 (01) 00841506137957	Screw, Non-Locking, 2.5mm x 26mm, Ti VAN-25260-TS (01)00841506138046	 (01) 00841506138046

Screw, Non-Locking, 2.5mm x 13mm, Ti VAN-25130-TS (01)00841506137964  <small>(01) 00841506137964</small>	Screw, Non-Locking, 2.5mm x 28mm, Ti VAN-25280-TS (01)00841506138053  <small>(01) 00841506138053</small>
Screw, Non-Locking, 2.5mm x 14mm, Ti VAN-25140-TS (01)00841506137971  <small>(01) 00841506137971</small>	Screw, Non-Locking, 2.5mm x 30mm, Ti VAN-25300-TS (01)00841506138060  <small>(01) 00841506138060</small>
Smooth Peg, Variable Angle Locking 1.8mm x XXmm, CoCr	
Smooth Peg, Variable Angle Locking, 1.8mm x 6mm, CoCr SPLS-18060-CC (01)00841506140360  <small>(01) 00841506140360</small>	Smooth Peg, Variable Angle Locking, 1.8mm x 13mm, CoCr SPLS-18130-CC (01)00841506140438  <small>(01) 00841506140438</small>
Smooth Peg, Variable Angle Locking, 1.8mm x 7mm, CoCr SPLS-18070-CC (01)00841506140377  <small>(01) 00841506140377</small>	Smooth Peg, Variable Angle Locking, 1.8mm x 14mm, CoCr SPLS-18140-CC (01)00841506140445  <small>(01) 00841506140445</small>
Smooth Peg, Variable Angle Locking, 1.8mm x 8mm, CoCr SPLS-18080-CC (01)00841506140384  <small>(01) 00841506140384</small>	Smooth Peg, Variable Angle Locking, 1.8mm x 15mm, CoCr SPLS-18150-CC (01)00841506140452  <small>(01) 00841506140452</small>
Smooth Peg, Variable Angle Locking, 1.8mm x 9mm, CoCr SPLS-18090-CC (01)00841506140391  <small>(01) 00841506140391</small>	Smooth Peg, Variable Angle Locking, 1.8mm x 16mm, CoCr SPLS-18160-CC (01)00841506140469  <small>(01) 00841506140469</small>
Smooth Peg, Variable Angle Locking, 1.8mm x 10mm, CoCr SPLS-18100-CC (01)00841506140407  <small>(01) 00841506140407</small>	Smooth Peg, Variable Angle Locking, 1.8mm x 18mm, CoCr SPLS-18180-CC (01)00841506140476  <small>(01) 00841506140476</small>
Smooth Peg, Variable Angle Locking, 1.8mm x 11mm, CoCr SPLS-18110-CC (01)00841506140414  <small>(01) 00841506140414</small>	Smooth Peg, Variable Angle Locking, 1.8mm x 20mm, CoCr SPLS-18200-CC (01)00841506140483  <small>(01) 00841506140483</small>
Smooth Peg, Variable Angle Locking, 1.8mm x 12mm, CoCr SPLS-18120-CC (01)00841506140421  <small>(01) 00841506140421</small>	Smooth Peg, Variable Angle Locking, 1.8mm x 22mm, CoCr SPLS-18220-CC (01)00841506140490  <small>(01) 00841506140490</small>